



Euronda®

It's all about protection

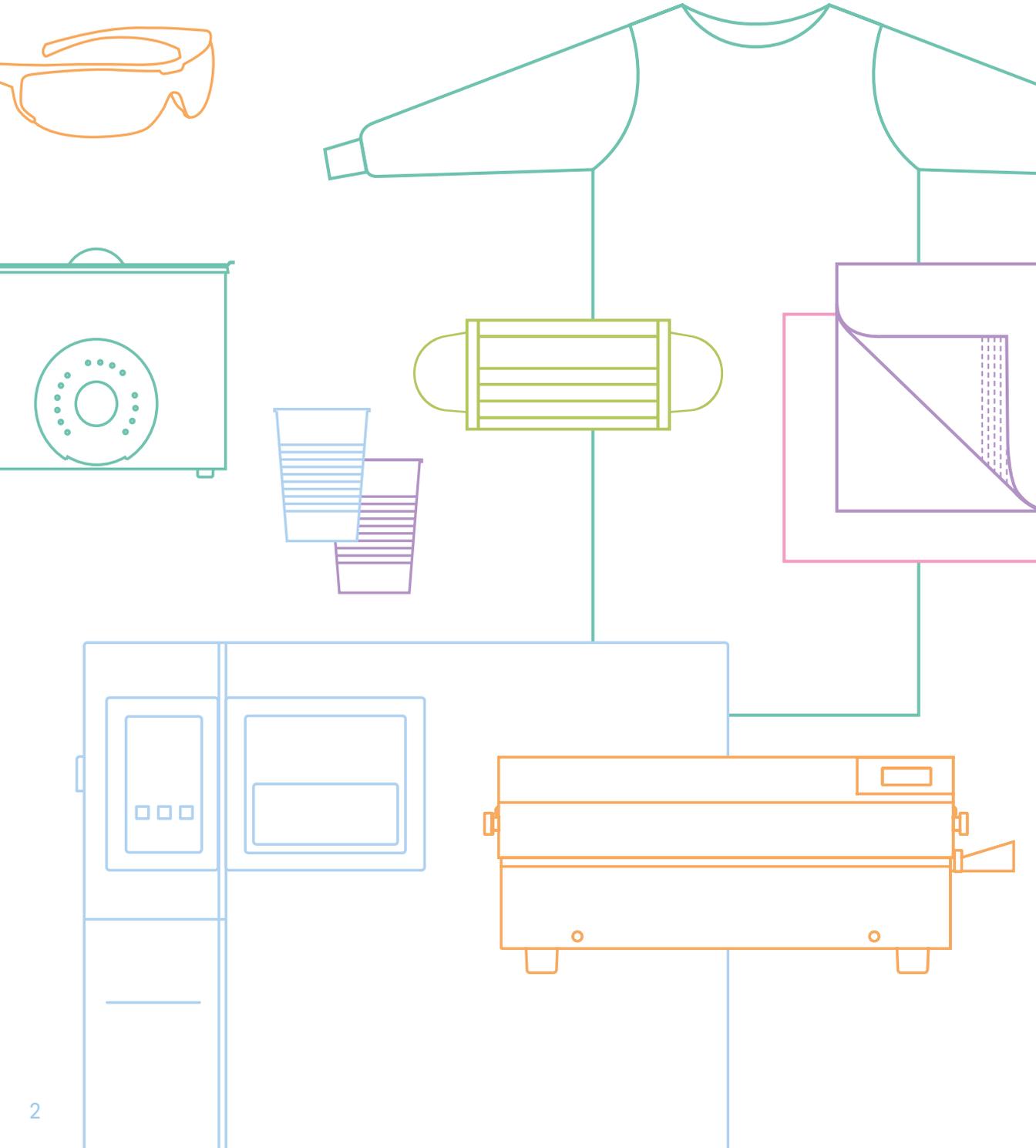

Index

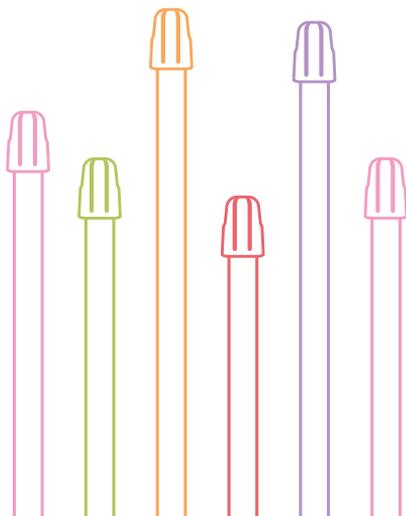

Monoart® News

- 4 **Protection 4 Face Masks**
Sensitive
- 5 **Barrier Film**
- 6 **Kit Monoart 100%**
The original!
- 7 **Plastic Cups**
New material
- 8 **Burgundy**
New color
- 9 **Blue Lagoon**
New color

Pro System News

- 10 **Eurosonic® 4D**
A new dimension
to ultrasonic washing
- 11 **Eurosonic® 3D**
The balanced choice
for ultrasonic washing
- 11 **Eurosonic® Energy**
Safe and accessible
ultrasonic washing
- 12 **Euromatic®**
Simplifies, speeds up
and keeps track of work
- 13 **Euroseal® Valida**
The thermosealer that
validates the process
- 13 **Euroseal® Infinity**
A sealing machine
is forever
- 14 **EXL**
Extra Large autoclave

Alle® News

- 16 **Doctor Set**
- 16 **Drape with U cut
without adhesive
and adhesive point**
- 17 **All-purpose bag
for radiographic systems**
- 17 **Medical adhesive**
- 18 **Mechanical irrigation lines**

Monoart[®] Protection 4 Face Mask Sensitive



The new Protection 4 Sensitive Face Masks are ideal for anyone seeking full operational protection and the safety of hypo-allergenic materials as well.

Made from hypo-allergenic cellulose, devoid of latex, glass fibre and chemical substances, they protect sensitive skins and avoid the risk of irritation, dermatitis or allergies.

They guarantee optimal protection against dust, powder, bacteria, aerosol spray and fluids, but they also feature a brand-new fitting system as well as a hydrophilic layer against the skin which absorbs sweat and moisture, offering improved comfort levels and Perfect fit.

Item code: 217141



Classification

Medical Device Class I
EN 14683: Type II
ASTM F2100: Level IIR



Color

Light Blue



Packaging

Box of 50 masks with an
“open/close” system
Carton of 20 boxes



CE Type Ref

Monoart[®] FACE MASK
PROTECTION 4
SENSITIVE



Perfect fit: the new “darted” system that makes for a better closer fit to the face.



Nose bridge extended to 130 mm for greater stability, externally welded elasticated ties so it lies flatter against the cheeks.



New package with an “open/close” system offers maximum practicality of use and preserves the integrity of the masks.

Monoart® Barrier Film



Disposable film for the protection of areas and surfaces, especially those difficult to reach, clean or disinfect. The light adhesive of the barrier film adheres easily to surfaces without leaving any residue.

The barrier film has a non-adhesive edge on both sides, to simplify placement and removal. The use of the barrier film reduces the risk of cross-contamination caused by aerosol and spraying. The pack contains 1 roll of 1200 tear-off 100x150 mm sheets with dispenser.

Item code: 269001



Classification

DM Class I



Dimension

100x150 mm



Color

Blue



Packaging

Single pack



CE Type Ref

MONOART® BARRIER
FILM



Kit Monoart® 100%

The original!



Available from
May 2019



Many colours, 100 items and 100% original. All the essence of Monoart®. All the original appeal of a unique range. All the quality of a 100% Italian product. All the safety and protection professionals and patients deserve. All the attention to detail typical of the Euronda brand.

The idea of this kit is to allow the dentist to have more colours, or to try out new colours, with a small investment. The best solution to try out 7 Monoart® colours and 5 products for yourself.

Practical Kit* containing:

- 100 Monoart® Towel Up!
- 100 Monoart® Saliva ejectors EM15
- 100 Monoart® Plastic cups
- 100 Monoart® Face masks PTC3
- 80 Monoart® Aprons PG30

** For the item code of the Kit Monoart® 100%, please refer to the Euronda Executive Assistant.*



Colors

Lime, green, blue, light blue, lilac, orange, blue lagoon**
** Available from July 2019



Box dimensions

31x31x24 cm



Weight

3,2 kg



Packaging

1 carton containing 1 Kit
Pallet with 64 boxes
120x80x200cm - 218 kg



Colourful

The pleasure of changing colour, the freedom to mix and the simplicity of trying something new.



Smart

Obtain maximum variety and quality at an affordable price and without inventory investments.



Compact

It offers maximum product rotation for the dentist, occupying minimum space, making the most of every centimetre.

Monoart[®] Plastic Cups New material



Monoart[®] disposable cups have been made with a new material and the polystyrene has been replaced with polypropylene.

All the new polypropylene cups will bear an embossed company logo, testifying to their high quality and safety standards Euronda. In addition, each box has a barcode label for a better storage facility.

Available in numerous colours and with floral pattern to customize your clinic.

Item codes: fuchsia **21410022**, floral pink **21410023**, pink **21410016**, red **21410017**, orange **21410018**, yellow **21410014**, lime **21410020**, floral lime **21410025**, green **21410013**, blue lagoon **21410047**, floral light blue **21410024**, light blue **21410012**, lilac **21410021**, floral lilac **21410026**, blue **21410015**, black **21410019**



Capacity

200 cc



Colors

Fuchsia, floral pink, pink, red, orange, yellow, lime, floral lime, green, blue lagoon, floral light blue, light blue, lilac, floral lilac, blue, black



Packaging

Bag of 100 pieces
Carton of 30 bags



Ergonomic

With a rounded edge and a secure grip.



Capacious

Available in 166 cc (white only) and 200 cc (multicolour) versions.



Resistant

Made with a new material: polypropylene instead of polystyrene.

Burgundy

New color



Intense and passionate, Burgundy is an original and sophisticated colour that playfully blends mystery and seduction. Matched up with the right colours, it creates truly elegant and refined combinations that will grab anyone's attention. It's a warm colour that creates an intimate and intense work environment, perfect for dynamic characters committed to change.

The new Burgundy is available for the following products:

Monoart® Face Mask PTC3 Floral (item code: 217126)

Monoart® Aspirator tips EM19 (item code: 22930009)

Monoart® Aspirator tips EM21 (item code: 22931009)

Monoart® Towel Up (item code: 21810414)

Monoart® Apron PG30 - 80 aprons (item code:


22010373) - 60 aprons (item code: **22010379**)

Monoart® Saliva ejectors EM15 - 15 cm fixed

(item code: **22810645**) - 15 cm removable (item

code: **22810644**) - 12,5 cm removable (item code:

22810646)

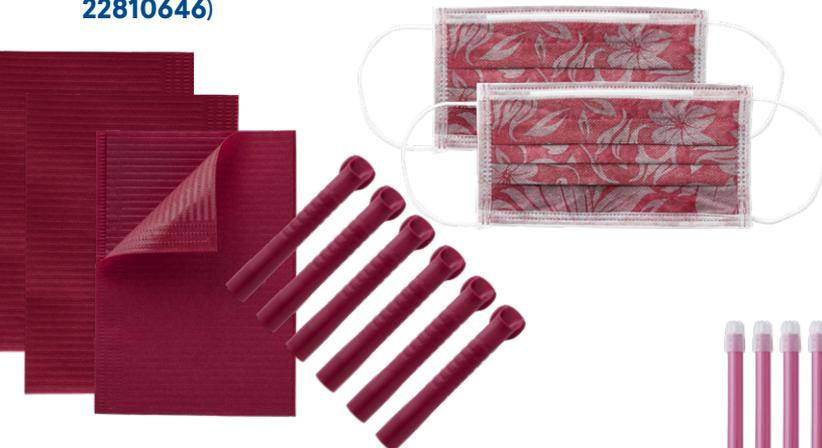

Blue Lagoon

New color

Elegant and sophisticated, Blue Lagoon combines two fluid shades of blue and green, embellishing them with an intense, material touch. Serenity blends seamlessly with freshness and harmony to create a beautiful nuance that pairs perfectly with light shades like white, as well as darker tones such as bordeaux and black, for a distinctive, refined effect.

The new Blue Lagoon is available for the following products:


Monoart® Plastic cups (item code: **21410047**)

Monoart® Towel Up (item code: **21810491**)

Monoart® Apron PG30 - 80 aprons (item code: **22010406**) - 60 aprons (item code: **22010407**)

Monoart® Saliva ejectors EM15 - 15 cm fixed (item code: **22810828**) - 15 cm removable (item code: **22810826**) - 12,5 cm fixed (item code: **22810827**) - 12,5 cm removable (item code: **22810825**)

Monoart® Face Mask PTC3 (available from July 2019)

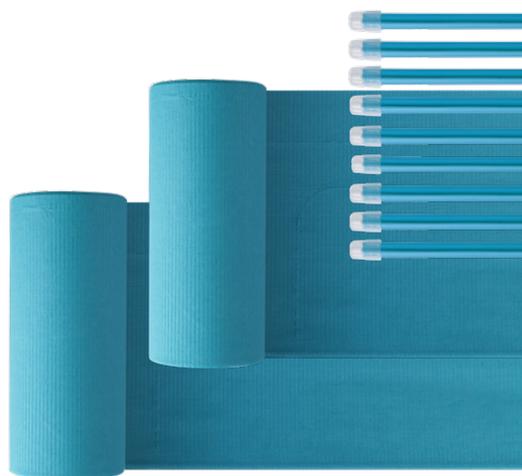

Eurosonic[®] 4D

A new dimension to ultrasonic washing

Eurosonic[®] 4D is the most advanced digital ultrasonic tank in the Euronda Pro System range. It has a LED backlit controller and it is intuitive and user-friendly thanks to the light signals associated with different cycles, times and temperatures.

There are two operational modes, allowing users to choose one of the 6 preset cycles or run a free cycle at a temperature between 30°C and 60°C, with a duration between 10 and 40 minutes. There are 4 different operating stages in the 6 preset cycles: Degas, Boost, Normal and Sweep.

Eurosonic[®] 4D has a drain on the back and it comes with a tap and pipe, as well as side handles for easier movement.

Item code: 113152



Tank capacity

3.8 L



Dimensions (wxhxd)

335x260x205 mm



Weight

5.8 Kg
(with basket and cover)



CE Type Ref

Eurosonic[®]: 4D



Eurosonic® 3D

The balanced choice for ultrasonic washing

A digital ultrasonic tank featuring a constant operating temperature of 60°C and time adjustable to between 0 and 30 minutes, through the backlit controller. Eurosonic® 3D is equipped with the Sweep Mode Technology for homogeneous distribution of the ultrasonic waves and more uniform cleaning. It also features the visual cycle countdown, thanks to the LED lights surrounding the front controller. The side handles facilitate the tank shifting operations.

Item code: 113102



Tank capacity

2.7 L



Dimensions (wxhxd)

270x240x195 mm



Weight

3.7 Kg
(with basket and cover)



CE Type Ref

Eurosonic®: 3D

Eurosonic® Energy

Safe and accessible ultrasonic washing

The new Eurosonic® Energy tank is equipped with Sweep Mode technology which distributes ultrasonic waves uniformly, encouraging thorough and even cleaning. Sturdy and affordable, with a steel body and anti-deformation non-drip cover, it has side handles to facilitate emptying and a new, more effective and intuitive dashboard. The cover Eurosonic® Energy is sound absorbing, as it is made of a special type of plastic, which attenuates noise produced by the ultrasonic waves.

Item code: 113025



Tank capacity

2.7 L



Dimensions (wxhxd)

270x240x220 mm



Weight

3.7 Kg
(with basket and cover)



CE Type Ref

Eurosonic®: Energy

Euromatic®

Simplifies, speeds up and keeps track of work

With a completely revamped design and new functions, Euromatic® is a rotary thermosealing machine with continuous cycle technology. It has a sealing speed of 8 m/min, which makes it ideal for clinics with a large number of pouches to seal. Data selected by the operator are printed directly onto the pouch by the built-in printer which features a new system for optimizing print cartridges. Not only does it ensure outstanding results in record time, it allows users to validate the process and keep track of all sealing data, also saving them on a USB device or SD card.

Item code: 108215



Sealing dimensions

12.5 mm



Dimensions (l x h x d)

510x215x242 mm



Weight

14 kg



Power supply

200-240 V | 50/60 Hz
Power consumption 600 W



**Available from
September 2019**



Euroseal® Valida

The thermosealer that validates the process

Euroseal® Valida comes with all of the best features of a manual sealing machine: it has an innovative and ergonomic design, it is sturdy, safe, and versatile thanks to its patented roll holder. Moreover, Euroseal® Valida can validate the sealing cycles and save data on a USB key or SD card. Sealing data can be transferred directly to a pc using Ethernet and WiFi connections. The 4" touch-screen display has a new user-friendly interface thanks to the upgraded software. Inserted directly on the support surface, users can quickly and intuitively browse through the menu and select one of the preset cycles or a free cycle.

Item code: 108301



Sealing dimensions

12x310 mm



Dimensions (lxhxd)

500x315x490 mm



Weight

8.3 kg



Power supply

200-240 V | 50/60 Hz
Power consumption 100 W




Available from
June 2019

Euroseal® Infinity

A sealing machine is forever

An electronic thermosealing machine for continuous cycle processing at a constant temperature: Euroseal® Infinity is the natural successor to Euroseal® 2001 Plus, from which it has taken and built upon its key elements. Robust, reliable and easy to use, the new design is both eye-catching and functional. The new built-in cutter and sealing lever positioned at the front not only improve ergonomics and safety but also make packaging and sealing operations easier. This design makes the machine more compact so it can be comfortably placed on a table top or wall-mounted in a vertical position.

Item code: 108145

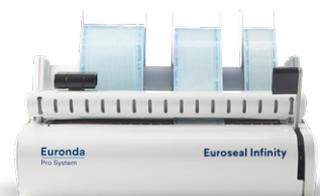

Sealing dimensions

12x310 mm



Dimensions (lxhxd)

430x215x382 mm



Weight

7.1 kg



Power supply

200-240 V | 50/60 Hz
Power consumption 100 W



Available from
October 2019



Capacity

With its 29 litre chamber, EXL is the biggest-capacity Euronda autoclave available, whilst ensuring optimum sterilisation results.



Traceability

It provides an expandable system of traceability thanks to standard connectivity and integrable Print Sets for the printing of adhesive labels.



Reliable

Based on the most highly developed Euronda autoclaves, EXL combines their strengths and best features in terms of duration, maintenance and design build.



Connectivity

The integrated SD Card and the Ethernet or WiFi connection will ensure that sterilisation data is always available, transferring it from the autoclave to a computer or other devices.



Eco-friendly

Thanks to its instant steam generator and the Pro System technology, sterilisation settings are adapted to the effective load, thus optimising water and electricity consumption.



Functional

From its user-friendly interface, to the most innovative systems which guarantee the autoclave's performance over time: EXL improves and simplifies sterilisation and the management of the sterilisation process.



EXL

Extra Large autoclave



The EXL was developed to meet the needs of clinics with high patient volumes, thanks to its greater load capacity, improved performance and low consumption. The large 29 litres chamber can contain up to five 55 cm trays, reducing the number of cycles necessary. The instant steam generator and the advanced drying cycle ensure optimum performance, with average water consumption of just 0.7 litres per cycle. The E-Touch colour display allows for a quick, intuitive use of the menu, and all operations can be smoothly managed thanks to the cycle programming functions (E-Timer). Traceability is guaranteed thanks to the SD Card and the possibility, on request, for two different types of print sets to be installed. The connectivity (Ethernet) and the software provided (E-Memory, E-Data) ensure results are available on a variety of supports and devices, so that all data can be shared and managed safely and in real time. EXL is designed to maintain top performance at all times, thanks to the integrated conductivity meter, the air-water separator and the air intake openings on the body. The magnetic top lid and the removable side panels make cleaning the tank easier, quicker and more effective.

Item code: 115760



Classification

Medical Device Class IIb
The Class B Autoclave fully conforms to EN 13060 standard requirements. Stainless steel moulded chamber certified to PED 2014/68/EC



Chamber capacity

29 L



Weight

62 Kg (machine empty with basket and trays)



External dimensions (wxhxd)

460x455x855 mm



Power consumption

3200 W 14 A



Water consumption

700 ml




CE Type Ref

EXL



Available from
August 2019

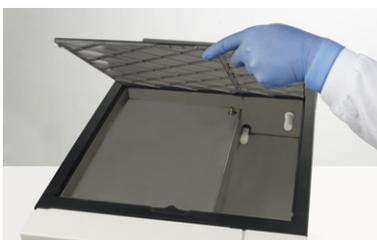

Easy access to the clean water tank



EXL is prepared for 5 trays or 3 shelves for cassettes



E-Touch Display

Doctor Set



Set designed specially with the dentist in mind, for brief and medium-length procedures. New composition with gowns in SMS microfibre, guaranteeing greater comfort for professionals.



Sterility

Sterility guaranteed for 5 years from the date of production

Set Code	Content	Size	Material	Pcs
270126	Microfiber gown with cotton cuffs	L	Non-woven fabric	1
	Surgical cap with adjustable ties	One size	Non-woven fabric	1
	Monoart® face mask PTC3	One size	Non-woven fabric	1
	Waterproof shoe covers	One size	CPE	2

Drape with U cut without adhesive and adhesive point



The 'U' shape makes it easy to place on the bottom lip or chin. The adhesive part allows you to fix the drape in place, which prevents it from moving. Industry experts maintain that its 115-cm length is perfect for covering patients. Packed in a sterile medical blister pack. The double adhesive label with barcode means that the product is traceable.



Classification

DM class I s



Material:

Non-woven fabric + polyethylene



Color

Light Blue



Sterility

Sterility guaranteed for 5 years from the date of production

Code	270230
Pcs	40
Sizes (wxl)	75x115 cm
Opening sizes (wxl)	15x25 cm
Adhesive point dimension	6x6 cm

All-purpose bag for radiographic systems



This all-purpose bag is sterile, see-through and allows you to protect radiographic, photographic and microscopic systems, leaving your surgical field totally sterile.



Classification

DM class I s



Material

PLT



Dimension

Ø 75 cm



Sterility

Sterility guaranteed for 5 years from the date of production

Code	Pcs
270316	25

Medical adhesive



Medical adhesive for securing protective sheaths, tubes and drapes, leaving your surgical field safe and sterile.



Classification

DM class I s



Material

Non-woven fabric



Dimension

10x5 cm



Sterility

Sterility guaranteed for 5 years from the date of production

Code	Pcs
270315	50 blisters of 2 pcs

Mechanical irrigation lines

An appropriate, well gauged jet of liquid coolant is able to bring a dental drill “to the ideal temperature”. This can only be done with an irrigation system that has been specifically created for the different types of physiodispenser available on the market. The soft silicone tips of Alle® irrigation lines can be adapted to the different diameters of wand cooling needles.



Classification

DM class II a



Material

Silicone, PVC, Phthalate-free ABS, Latex-free



Sterility

Sterility guaranteed for 5 years from the date of production

*All brands mentioned are registered and owned by the producers.

Details

Code	Type	Models*	Pcs
270601 NEW 2019	1/2Y	Suni Satelec 2000/3000/Suni Max Surgi set basic, Surgi set pro 500, Dental surgery unit, Dental Unit DSC Electronic, Sweden&Martina XO Osseo, Sky unit	10
270602	1/2Y	W&H Elcomed 100/200 / Astratech Elcomed 100	10
270604	1	Nouvag TCM3000, MD8000, MD7000, MD10, MD10S, MD11, MD20, MD30	10
270606	1/2Y	De Giorgi Steril Intraplant, Intrasurgery 2.0, Intramatic 2, Intramatic Plus	10
270607	1/2Y	Kavo Intrasept 905 / Mectron Piezosurgery II, Piezosurgery 3, Piezosurgery Touch	10
270609	1/2Y	NSK Surgic XT, Surgic XT Plus, Osseocision®, VarioSurg / Aseptico Aeu 6000-70V, 7000E-70V, 7000L-70V, 17Bv2, 26OS, 707Av2 / Bonart Piezo Surgical System	10
270610	1/2Y T	W&H Implantmed, Elcomed SA-310, Piezomed / Nobel Biocare OsseoSet™	10
270613	1/2Y	Anthogyr Implanteo, Implanteo Led	10
270616	1/2Y	W&H Elcomed SA-200C / Frios Unit E for Dentsplay	10
270617	1/2Y	KaVo Intrasurg 300/500, Intrasurg 300 PLUS, Intrasurg 1000, Intrasurg, 1000 Air	10
270618	1/2Y	Bien Air Chiropro 980 / BioSAFin Easy Surgery, MyTUTOR / EMS Piezon, Master Surgery® / Esacrom Surgysonic I Surgysonic II, Surgysonic, Moto	10
270621	1	Kavo MASTERsurg, EXPERTsurg	10
270622	1	Nouvag Physiodispenser 2000	10

Code	Models*	Pcs
270626 NEW 2019	Extension line for irrigation set with luer lock connector	10

Legend

1 - Single irrigation

1/2 - Internal and external irrigation without a Y-connector

1/2Y - Internal and external irrigation with a Y-connector

T - Tap to connect a second irrigation line to a second handpiece or to connect an irrigator

Euronda Spa

Via Chizzalunga 1
36066 Sandrigo
Vicenza, Italy
t (+39) 0444 656111
f (+39) 0444 656199
m info@euronda.com
www.euronda.com

Deutschland

Euronda Deutschland GMBH
Am Landwehrbach 5
48341 Altenberge, Deutschland
t (+49) 2505 9389 0
f (+49) 2505 9389 29
m info@euronda.de
www.euronda.de

France

Euronda France Eurl
ZAC les Vallées, Avenue de Bruxelles
60110 Amblainville, France
t + 33 (0) 3 44 06 69 70
f + 33 (0) 3 44 02 03 89
m info@euronda.fr
www.euronda.fr

España

Euronda Spa
Identyd, Av. Cabezo de Torres 4
30007 Murcia, España
t 900102034
m info@euronda.es
www.euronda.es

Russia

OOO Finstar
5/2a Varshavskaya str.
196128 St. Petersburg, Russia
t + 7 (812) 635 88 94
m info@finstartech.ru
www.euronda.com.ru

Euronda